




APRIL 2019

MYOGLOBIN ASSAY

Our Myoglobin test is based on the principle of latex particle-enhanced immunoturbidimetric measurement. Latex particles are coated with human myoglobin antibodies, forming complexes with the antigen present in the sample.

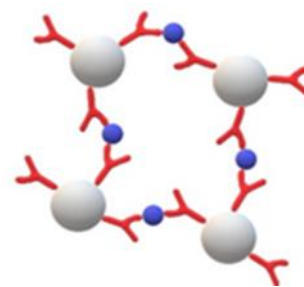

General information:

Sample material	Serum, Plasma (EDTA, heparin, citrate)
Shelf life	18 months from date of production
Method	Immunoturbidimetric
Reaction	Fixed time
Wavelength	580 nm
Linearity	5 - 600 ng/mL
Limit of Detection	5 ng/mL

Applications for automated clinical chemistry analyzers are available upon request!

Order information for the Myoglobin assay:

REF	Content	
909240	2x 27 mL Buffer + 2x 9 mL Latex Reagent	Myoglobin
909260	4x 1 mL	Myoglobin Calibrator 4 Level Series
D18270	2x 1 mL	Myoglobin Control Set (2 Levels)

**Your
DIALAB-Team**

